

DARCO Body Armor Heel Reliever

Effective pressure relief in the heel area, especially for supine patients.



The Body Armor® Heel Reliever redistributes heel pressure to the gastrosoleus complex in the calf. The heel area is completely free of all contact pressure. The heel “floats” freely in a pressure-free zone.

Forefoot strap and forefoot pad
 > keeps the foot safely in proper position
 > for pressure redistribution
 > forefoot pad prevents skin contact with forefoot strap

Foam with bevelled edges and corners
 > prevents pressure at the edges
 > facilitates the overlapping closure

Welded seams
 > for less irritation

Hook and loop straps
 > for easy handling
 > secures the foot in a comfortable position
 > with sewn stitched edges

Ventilation holes
 > promotes air circulation
 > helps to regulate the skin's environment

Product label
 > can be written on
 > comes with a QR code containing instructions for use

Extra tall foot section
 > accommodates larger feet
 > ensures an orthograde position
 > prevents foot drop
 > prevents foot rotation

Long calf pad (30 cm) fixed to the inside
 > transfers the load from the heel
 > avoids hyperextension of the knee

Heel opening
 > for visual control of correct heel suspension

Smooth, friction-free exterior surface
 > allows the foot to move freely in bed without affecting the suspended position of the foot

DARCO Body Armor Heel Reliever


- For the prevention and healing of heel pressure ulcers as well as diabetic foot syndrome and decubitus
- Braden scale scores of fourteen or less
- Diabetes with neurological sensory deficit
- Pedal pulse not palpable
- Paraplegia
- Prolonged immobility (> hours) or coma

- Product includes:
- Body Armor Heel Reliever
 - Forefoot Pad
 - Extra Customisable Foam Pad
 - Positioning Wedge
 - Laundry Bag


■ ONE SIZE FITS ALL

DARBAHR-F Flat Foam Interior
DARBAHR-N Napped Foam Interior

AVAILABLE IN TWO INNER LINER FOAM OPTIONS:

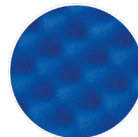

Napped foam
 results in optimum air circulation and helps to ensure a pleasant cooling effect.

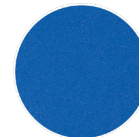

Flat foam
 developed especially for sensitive skin. It is particularly suitable for Diabetes patients. This helps to prevent the formation of oedema.

briggatemedicalcompany

Supporting Podiatrists with integrity and trust since 1987

(03) 8586 7800 (Melb. Metro)

1800 33 4142 (Free Call)

briggate@briggate.com.au

www.briggatemedical.com