

Ezy-As Stocking/Garment Applicator



Ezy-As Stocking/Garment Applicator

Designed to assist users with application of compression stockings, tubular bandages or garments. The compression stocking or tubular bandage is turned inside out over the applicator evenly distributing elastic forces of the garment onto the smooth outer surface of the device. Once the garment is correctly positioned on the applicator the user then transfers the compression garment onto the limb.

Features:

- Works with both open and closed toe stockings
- Fits compression garments to upper and lower limbs
- Significantly places less stress and pain on the injury and / or limb
- Robust, durable and easy to clean
- Compact, lightweight and travel friendly
- Australian Made and designed
- Available in three colour coded sizes for easy identification
- Attachable handles (sold separately) assist those with difficulty bending over


EZYAS-S	SMALL	Orange	Calf Measurement: 26cm-34cm
EZYAS-M	MEDIUM	Yellow	Calf Measurement: 32cm-40cm
EZYAS-L	LARGE	Blue	Calf Measurement: 38cm-47cm

Ezy-As Handles - sold separately

- Attachable handles assist those with difficulty bending over
- Allows patient to pull the device up over the area with ease

EZYAS-HAND	Attachable Handles	Fits all sizes
-------------------	--------------------	----------------

Log-in to www.briggatemedical.com for pricing

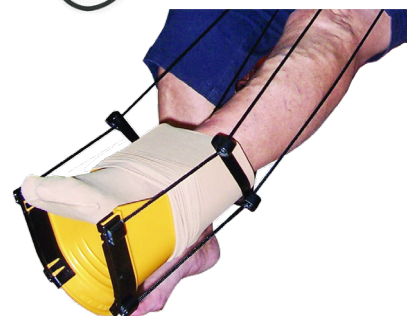

briggatemedicalcompany

Supporting Podiatrists with integrity and trust since 1987

(03) 8586 7800 (Melb. Metro)

1800 33 4142 (Free Call)

briggate@briggate.com.au

www.briggatemedical.com