



Attaches to the side  
of the Melag Euroklav  
23S+ or 23VS+.



## Melag Meladem 40 Water Management System

- Attaches to the right-hand-side of your Melag 23S, 23VS or 29VS autoclave
- Direct feed of demineralised water
- Allows automatic drainage of used/grey water
- Filter cartridges (on average) treat 300 Litres of tap water  
(this can vary depending on your local water quality)

**MELDEM40** MELAdem 40 Water Treatment Unit

**MELDEMFLT** Replacement Filter Cartridges Pair

Log-in to [www.briggatemedical](http://www.briggatemedical) for pricing

MELAdem 40 Filter Cartridges need to be replaced after 300 litres has been produced (Approx. 1 litre is used per cycle).

### System includes:

- Meladem 40 Unit
- 2 x Pairs of cartridges
- 2 x Brass Fittings (for 20mm BSB Tap & rear of Autoclave)
- 1 x 2m Drainage Hose with Connection

### Specifications:

- Dimensions: 32cmW x 35cmH x 15cmD
- Weight: 3.7kg (incl. Resin Filter Cartridges)
- Compatible with: Euroklav 23S+ / 23VS+ / 29VS+

Manufactured in  
Germany



### IMPORTANT INSTALLATION INFORMATION:


As this device is plumbed to mains water, a water conductivity test of your water supply is recommended to ensure compatibility. Failure to do so may result in damage to your Melag devices.

## Melag Melajet

- Optional extra that connects to the Meladem 40
- Water spray gun that assists in the cleaning and rinsing of instruments

**MELAJET** Melajet Spray Attachment for Meladem 40

Log-in to [www.briggatemedical](http://www.briggatemedical) for pricing

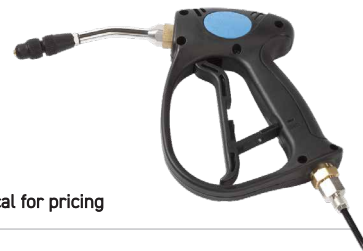

# briggatemedicalcompany

Supporting Podiatrists with integrity and trust since 1987