

Water Distiller

User Manual



IMPORTANT:

PLEASE READ THIS INSTRUCTION MANUAL BEFORE OPERATING THE WATER DISTILLER FOR THE FIRST TIME.


FAILURE TO DO SO MAY RESULT IN MACHINE MALFUNCTION, REDUCE THE LIFE-SPAN OF THE MACHINE AND/OR VOID WARRANTY.

TO PROLONG THE WORKING LIFE OF THIS UNIT, WE SUGGEST THAT AN EXTERNAL TIMING DEVICE (NOT INCLUDED) BE USED TO AVOID THIS DEVICE BOILING DRY DURING OPERATION.

BEFORE USE:

Unpack the Water Distiller and check all the components are present and undamaged. The box should contain:

- A. 1 x Water Distiller
- B. 1 x Water Storage Cannister
- C. 1 x Discharge Cap
- D. 1 x Cannister Cap
- E. 1 x Filter
- F. 1 x Mains Power Cord
- G. 1 x Jar of Chamber Cleaning Crystals

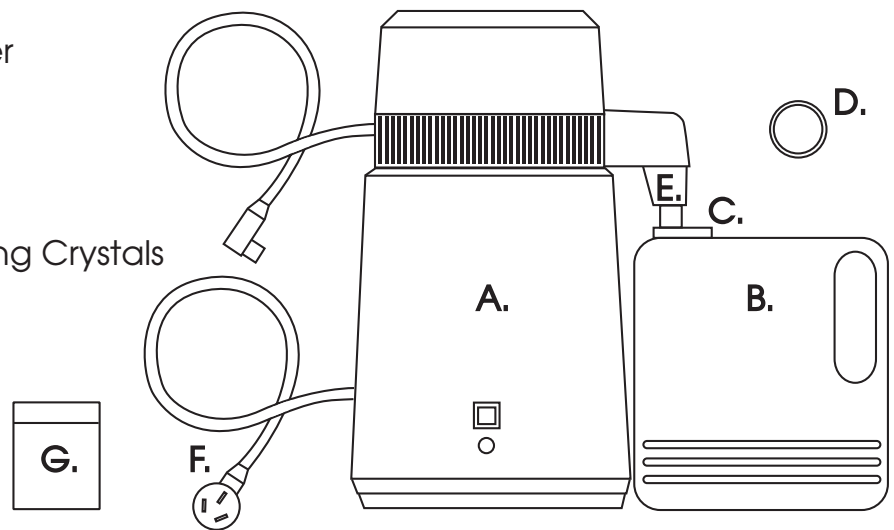

INITIAL SET-UP:

1. Remove Water Cannister from the inside of the Water Distiller.
2. Fill the Water Distiller Chamber with water to the highest point indicated (marked 'FULL' on the inside of chamber). You can reduce the time needed to distill water by pre-boiling the water to be distilled.
3. Place the Top Cover onto the machine's body whilst making sure the silicone seal is correctly in place around the entire rim. Ensure the unit is situated on a stable surface then connect the power cord from the Top Cover to the Main Unit Body.
4. Place the Discharge Cap on the Storage Cannister. Place the Cannister under the Filter. Connect the Mains Power Cord to power outlet.
5. Press the "Start Button". The fan will switch on automatically and the tank will begin to heat, thus starting the distilling process.
6. Distilled water will begin to drip into the Storage Cannister after approximately 30 minutes. It will take approximately 4 hours to distill 4 litres of water. The unit will automatically shut down when the cycle has finished.

PLEASE NOTE:

The unit must be allowed to cool down before starting another distilling cycle (this will help prevent injuries from scalding water and prolong the life-span of the Water Distiller).

continued next page

DO NOT USE THE FIRST CANNISTER OF DISTILLED WATER PRODUCED BY THE MACHINE. IT TAKES ONE FULL DISTILLING CYCLE FOR THE UNIT TO 'SELF-CLEAN'. ONCE THE UNIT HAS COOLED, RINSE AND WIPE OUT THE TANK. DISCARD WATER PRODUCED FROM FIRST CYCLE.

IMPORTANT NOTES:

- **Disconnect Mains Power Cord when unit is not in use. When Mains Power Cord is plugged in, always ensure there is water in the chamber.**
- Use the unit in a dry, well ventilated area.
- Keep the body of the unit dry - only fill with water using a jug.
- "Start Button" will begin distilling cycle and the unit will automatically turn itself off once cycle has finished. To stop the unit mid-cycle, remove Mains Power Cord from socket.
- Do not immerse any part of the unit in water.
- Do not remove the Top Cover whilst unit is mid-cycle. Wait until unit has cooled down before removing top.

CLEANING:

Rinse and wipe-out the tank chamber everyday using clean water. Calcification and scale will build up on the bottom of the tank and should be removed approximately every 2 weeks (depending on frequency of use). To clean the tank, follow the instructions below:

1. Pour enough hot water to cover scale in the tank.
2. Place 1-4 teaspoons of Cleaning Crystals into the tank depending on the amount of scale.
3. Leave the lid of the unit open and start up a cycle by pressing the "Start Button".
4. Let cycle run for 30 minutes, pour out water and wash the inside of the tank with clean water.
5. If scale is not completely removed - repeat the process.

FILTER REPLACEMENT:

The Water Distiller Filter should be replaced approximately every 2-3 months. To replace the filter, follow the instructions below:

1. Remove Top Cover and place upside-down.
2. Squeeze filter in towards Top Cover to release from water outlet.
3. Remove Filter from water outlet.
4. Check that new Filter is free of visible contaminants.
5. Insert new Filter using reverse of the above instructions.

briggatemedicalcompany

23-25 Lakewood Blvd
Braeside Vic 3195
Australia

Ph: (03) 8586 7800
Fax: (03) 8586 7888

briggate@briggate.com.au
www.briggatemedical.com